

Reducing compliance costs in the life sciences industry

Eurotherm®

Sterilization

Eurotherm™ provides a range of products, digital engineered solutions and services throughout the world. Our expertise in life sciences allows us to supply solutions that suit the scale of your manufacturing requirements, while helping to maximize efficiency, productivity, and ultimately your return on investment.

Our solutions support regulatory compliance and help create a safer world.

Compliance

We empower our customers to help maintain regulatory compliance and minimize audit costs by providing a Data Integrity layer with open IoT platforms that support the digital transformation to Pharma 4.0.

Safer world

Specialists in managing critical data and contextual metadata to efficiently manage the quality, safety and authenticity of manufactured goods.

We have application expertise in:

- Control and sequencing
- Recipe management
- Batch control and reporting
- Setpoint programming
- Graphics based on situational awareness concepts
- Alarm management
- FDA 21 CFR Part 11 and EudraLex Annex 11 compliance
 - Electronic signatures
 - Electronic records
 - Audit trail
- Data Integrity ALCOA+ concept
- ISPE GAMP® 5 guidelines
 - Good engineering practice (GEP)
 - Qualification practices
 - Risk based approach
 - Quality management
- System lifecycle support services

Sterilization overview

Sterilization is the act or process, physical or chemical, that destroys or eliminates all viable microbes including resistant bacterial spores from a fluid or a solid. Despite being stated as an absolute, the action of sterilization is usually stated in terms of probability of survival of a known quantity of a specific microorganism (F_0). There are many sterilization methods: chemical, dry heat, filtration, gas, steam, radiation, and vapor.

Heat sterilization



Steam heating: High-temperature moist steam is the most widely used and dependable sterilization process as it is relatively economical, nontoxic, fast acting on microorganisms, and offers adequate penetration. On certain materials there can be some downsides, including corrosion and overheating.

Dry heating: Used where the material could be damaged by high temperature or impenetrable by moist steam (e.g. powders, petroleum products, and sharp instruments). Dry heat is nontoxic, low cost, penetrates materials and is noncorrosive for metal and sharp instruments. Penetration and microbicidal action is generally slower than steam. Some machines are equipped with forced-air or mechanical convection systems that circulate air at a higher velocity and can reduce process time.

Gas and liquid chemical sterilization

Sterilization with gas is subject to risks of explosion and issues with gas toxicity. Common types are shown below:

- Ethylene oxide gas (ETO, EtO, EO) is typically used to treat items that are moisture or heat-sensitive. Performance qualification of the ETO sterilization process is complex because in addition to temperature, the humidity, vacuum/positive pressure, and gas concentration also require rigid control
- Chlorine dioxide gas has been shown to be effective at decontamination of large-scale facilities, rooms, isolators, processing vessels and tanks, and biological safety cabinets
- Vaporized hydrogen peroxide has been used as a chemical sterilant for many years. This method is currently being used in aseptic product filling isolators, sterility testing isolators, lyophilizers, and sanitization of clean rooms
- Vaporized peracetic acid systems are not currently cleared for use by the FDA

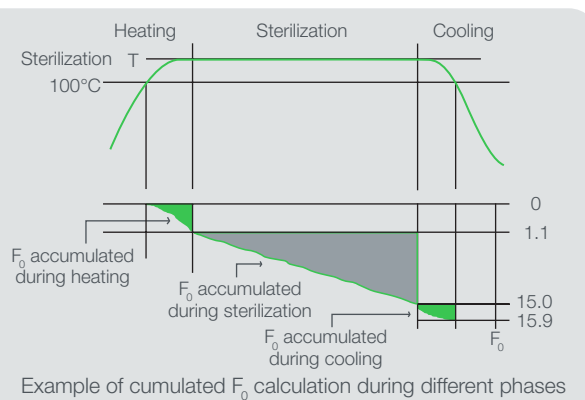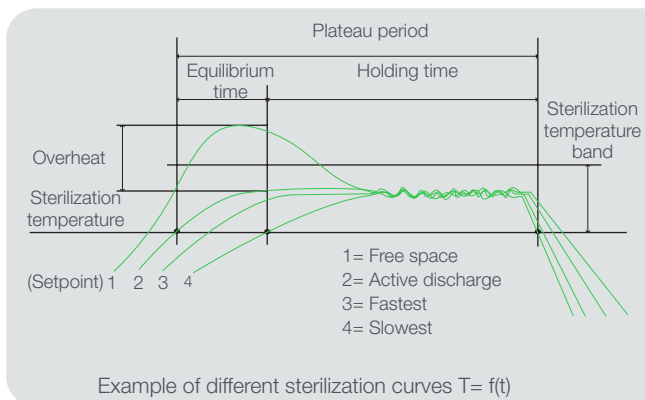

Steam process

The process is well understood and the attainment of sterilization conditions can normally be confirmed by simple physical measurements.

Sterilizers using high-temperature steam to process porous loads have a preferred sterilization temperature of 134°C¹. Aqueous “fluids sterilizers” operate at a preferred sterilization temperature of 121°C². Fluids sterilization is based on the definition of “probability of a non-sterile unit” (PNSU).

Dry heating processes

Dry-heat sterilizers are essentially electrical ovens. They operate around the sterilization temperature of 160°C³.

A filter and fan are used to maintain the chamber slightly above atmospheric pressure to ensure that the sterility of the product and the integrity of the cleanroom environment are not compromised.

Typical steps are: heat, sterilize, and cool.

UK NHS – Health Technical Memorandum 2010 Part 1 '2.11, '2.13, '2.20

Eurotherm solution:


- Distributed control system with configuration lock
- Fast acting control accuracy and process repeatability
- Functions for monitoring sterilization processes, including F0 calculation for parametric release
- Functions for calculating MKT (mean kinetic temperature) and steam flow
- Batch/recipe management
- Data management based on ALCOA+ principles
- Power control for electric heaters
- Local HMIs to full SCADA solutions
- High availability architecture (redundant solutions and 'Store and Forward' feature)
- Data analysis
- Historian and reporting server
- Reporting

FDA 21 CFR part 11 and EU EudraLex Annex 11 regulations compliancy

As an automation and information technology leader our solutions meet the Electronic Record and Electronic Signature requirements as defined by the US and EU regulatory bodies.

Data Integrity ALCOA+ guidelines

To make sound decisions, you need to trust your data. Key regulatory bodies (FDA, EMA, WHO) and some advisory bodies (PIC/S, ISPE) have agreed on the Data Integrity related ALCOA+ concept. ALCOA+ defines that data should be Attributable, Legible, Contemporaneous, Original, and Accurate + Complete, Consistent, Enduring, and Available. As an experienced solution supplier, we established in life science processes, Eurotherm is a main supporter of that vision and contributed to the definition and the revision of some of these guidelines.

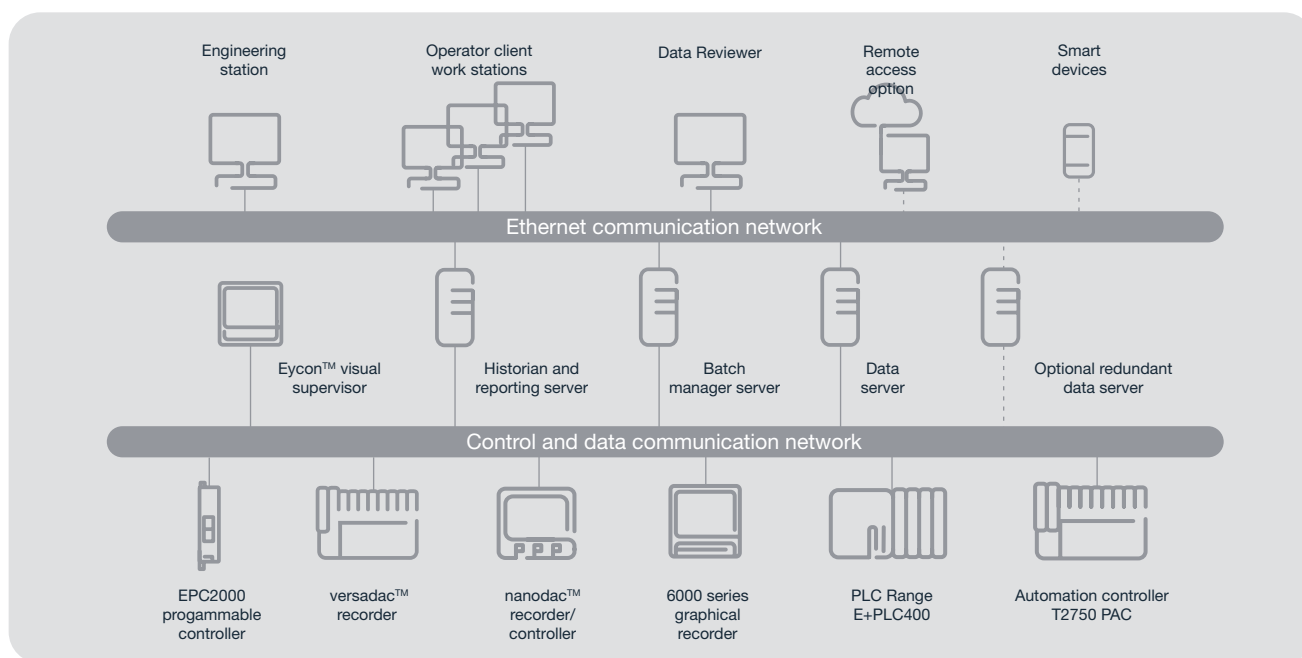

A 21st century risk-based approach

Business investments should be future-proof and audits hassle free. Eurotherm has developed and widely applied a set of good engineering practise (GEP) qualification documents based on ISPE GAMP 5 guidelines to assist with achieving these goals. Qualification documents can be maintained in electronic format. Industry is progressing from manual standard operating procedure (SOP) based manufacturing operations to a digitalized paperless quality systems approach based on FDA and ICH guidelines.

Quality by design

In a quality by design (QbD) approach, product quality is continuously monitored and controlled at the earliest stages, rather than waiting for a process to end. Pharmaceutical manufacturers need to focus on identifying, controlling, and validating process variables that could cause a non-compliant result. This is accomplished by managing the quality target product profile (QTPP), the critical quality attributes (CQA), and the critical process parameters (CPP). As defined by the process analytical technology (PAT) approach, Eurotherm can assist with measurement and performance analytics on CQAs and help manage CPP deviations, providing time-stamped evidence for correlation of parameter behaviors at their occurrence.

Pharma 4.0 ready technology

Eurotherm control and data recording solutions are IoT ready, providing a data integrity layer within open IoT platform system architectures and aiding the digital transformation to Pharma 4.0 technology.

Eurotherm US LLC

44621 Guilford Drive, Suite 100
20147 Ashburn, VA
USA

Phone: +1-703-724-7300

www.eurotherm.com



Contact your local sales representative



Document Number HA033627USA Issue 2

Watlow, Eurotherm, EurothermSuite, EFit, EPack, EPower, Eycon, Chessell, Mini8, nanodac, piccolo and versadac are trademarks and property of Watlow, its subsidiaries and affiliated companies. All other trademarks are the property of their respective owners.

©Watlow Electric Manufacturing Company. All rights reserved.

